


Aphelion Dev

IMAGE PROCESSING
AND IMAGE ANALYSIS SOFTWARE
TO FACE NEW CHALLENGES

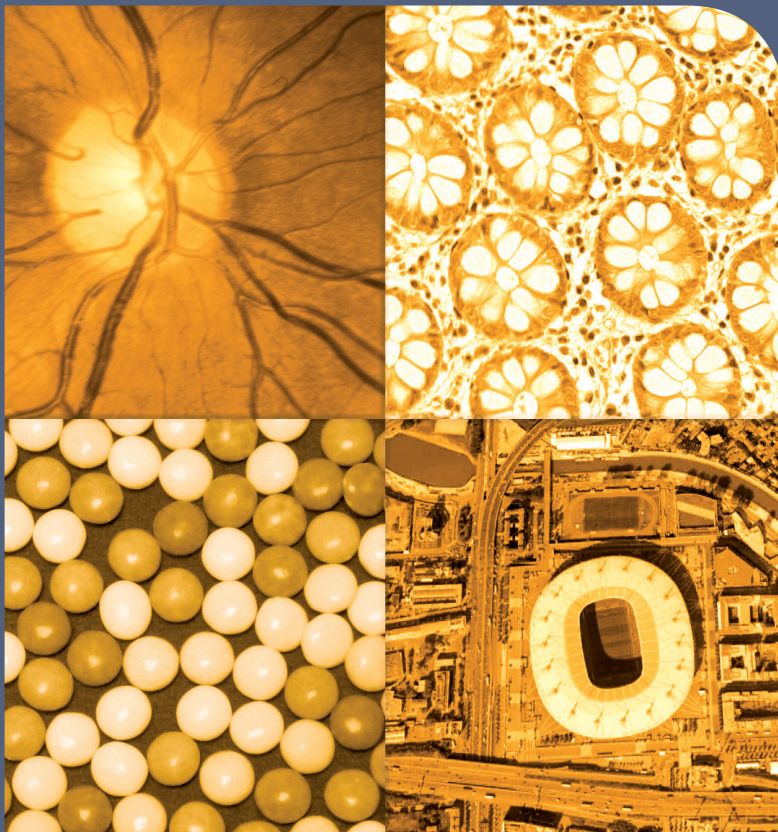

BIOLOGY / COSMETICS / GEOLOGY
INSPECTION / MATERIALS SCIENCE

METROLOGY
OBJECT RECOGNITION


PHARMACOLOGY / QUALITY CONTROL
REMOTE SENSING / ROBOTICS

SECURITY / TRACKING



Aphelion™ Dev to Quickly Prototype & Deploy Imaging Applications

We have been providing products to the image processing and analysis marketplace for more than 25 years. Our Aphelion software products and our customers benefit from that experience. While many competitor products have come and gone, Aphelion™ Dev remains the product of choice for those developing advanced image-based applications. Our Dev product has earned this position because we have continuously updated it to keep pace with imaging technology advances, and the evolution of Windows® environments and personal computer hardware technologies. At the same time we continued providing timely and professional customer support. Our keeping pace with rapidly changing imaging marketplace has substantially resulted from our partnering with international R&D labs and our successful track record developing advanced imaging applications.

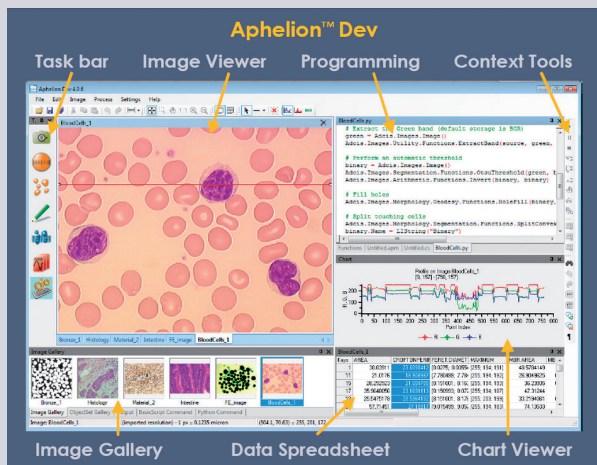

Aphelion Dev's graphical user interface utilizes the best features of current Windows® technology, such as .NET™ components and software architecture tools. Its context-based tools guide a new user from image acquisition, through image processing and analysis, to report generation. Yet, Dev provides the more experienced user the tools to develop advanced and automated imaging processes that can quickly become deployed applications.

Aphelion Dev users benefit from its openness, made possible by support for multiple macro languages that can call functions from Dev's libraries as well as from 3rd party software products. Dev's libraries are available as native code and .NET™ components. This capability enables OEMs and VARs to create stand-alone imaging applications with customized GUIs in a timely manner.

Keeping up with rapidly changing technology, Aphelion™ Dev supports monochrome, color, multiband, 2D, and 3D⁽¹⁾ images; a set of optional extensions for specific applications, and the most recent Windows® operating systems⁽²⁾. Aphelion Dev's design significantly enhances performance for both development work and for deployed systems. User productivity is also benefitted by an extensive Help library and video tutorials.




Settings

Settings are those parameters that configure the Aphelion Dev software for control of the image acquisition system being used, calibration of measurements for real-world units, and storage location of user data files, and many other parameters. Centralized control of these parameters facilitates the user making system changes.



Calibration Settings

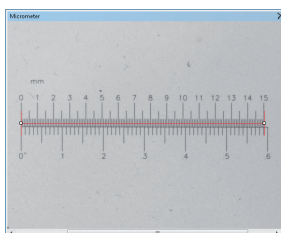

Aphelion Dev can be calibrated to report measurements accurately in real-world units. Using an object of known dimensions (e.g.,

micrometer) and Dev's on-screen measurement tools, a user can quickly create a calibration profile for a specific imager (e.g., camera) or a specific configuration of an imager. Dev does not limit the number of calibration profiles that can be saved, thereby enabling the generation of accurate, real-world measurements for any configuration of an imaging device that was calibrated with Aphelion Dev.



Measurement Settings

Measurement settings enable a user to select from a large library of measurements (e.g., length, area, texture, intensity statistics) that are to be used in an image analysis project for each predefined object types (e.g., line, circle, bitmap).


More than 50 different measurements can be computed per object, and there are no limit on the number of objects per image or on image size. If a calibration profile has been specified, all measurements are then provided in real-world units.



Image Acquisition

Aphelion Dev supports any acquisition devices compliant with the DirectShow® interface protocol (e.g., webcams and numerous cameras on the market) and the IIDC specification for an IEEE 1394 bus⁽³⁾ (i.e., FireWire® compatible). Optional interfaces broadens the range of compliant cameras. Once the proper driver for the digital camera is installed, a live image can be displayed in the Aphelion Dev interface on demand. The user has only to select the right camera configuration settings and adjust the exposure parameters, click the GUI's "Live" button to display a live image, focus the camera, apply a shading correction to get an even background, and then click the "Snap" button to capture a high quality image. Multiple camera settings can be saved and reloaded in the Dev.

Images previously captured and stored on media can also be loaded and analyzed with Dev. Thanks the 64-bit architecture support, very large images (e.g., from a scanner microscope or scanning electron microscope) can be processed, thus providing high throughput imaging.

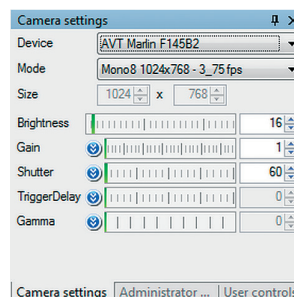




Object Extraction

Aphelion Dev provides two methods to extract objects of interest from an image. One method uses the Developer task to extract objects from complex images or to fully automate the extraction process. Dev includes a large set of extraction functions, from basic to very advanced. The second method uses Dev's Object Extraction task, a semi-automated process that guides the user through three steps: Color band selection, automatic or interactive thresholding, and object filtering to further refine the extraction.



Object Editing

The Object Editing task provides tools for modifying the set of objects that were extracted. Altering or removing existing objects and drawing new objects can be easily performed using mouse operations. An undo function is always available to go back to the previous step of the editing process. The combination of the Object Extraction and Editing tasks enables a user to achieve an optimal segmentation of an image regardless of its complexity.

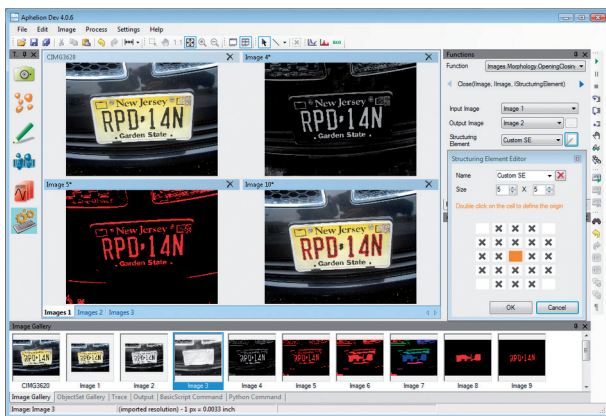

Measurements

Object Measurements

Aphelion Dev's processing of object measurements is very quick, and up to 50 measurements can be computed with just one mouse click. The computed measurements are displayed in a grid format for the objects automatically extracted and for those drawn with the Object Editing task. The GUI displays the grid, its associated image with its objects overlaid on the image, and an interactive chart displaying a histogram or scatter plot of object measurement values (i.e., values in a grid column). Any object can be highlighted simultaneously in both the image overlay, the chart, and the grid by clicking on the object in either the overlay or the object's row in the grid. An export feature saves grid data in Microsoft® Excel® or CSV format.

Interactive Measurements

Interactive measurement tools are used to draw virtually any geometric shapes in the objects overlay in order to quickly and directly measure object features (e.g., length, circumference, area). This is especially useful when there are only a few objects to be measured or when ad hoc measurements are desired. Measurements for the shapes drawn are displayed in a grid format as well as in the image overlay. As with the Object Measurements grid, the Interactive Measurements grid can be exported to

Excel® or CSV Text format. In addition, interactive histograms in a region of interest and intensity profiles are provided.



Developer

The Developer task is Aphelion Dev's primary tool for developing and deploying sophisticated, image-based applications. It provides hundreds of image and ObjectSet functions to enhance, restore, and filter images; extract objects such as blobs, chains, and lines; and transform, group, analyze, and classify objects. These can be applied to an entire image or to multiple regions of interest of any size and shape. Developer's functions are all available from drop-down menus. When a function is selected, Dev provides a dialog box that enables easy entry and full control of the function's parameters. Specialized dialog boxes facilitate entry of complex parameters (e.g., kernels). Once its parameters are entered, the function can be executed with a mouse click and results immediately displayed.

To create and automate complex image processing and analysis applications, Dev's functions can be combined using a macro language accessible without leaving the Dev GUI. Four macro languages are supported: BasicScript®, an easy to use scripting language, Python®, Visual C#®, and Visual Basic®. Debugging tools and a macro recorder⁽⁴⁾ are provided to aid in program development.

In addition to the comprehensive function libraries for image and object processing and measurement, programmers will appreciate the extensive set of utility functions included in Dev. These provide for image and data management, display controls, graphics handling, and support for many image formats. All Dev functions are fully documented with detailed descriptions of properties and methods, and with coding examples.

The Developer task has been designed for application development in virtually any image domain. It is especially well-suited for applications in the fields of biology, microscopy, materials science, quality control, and remote sensing.



Report Generation

The final step of a typical image analysis project is to generate a report of the results obtained from the analysis. In the past, this was often done by using multiple software products, requiring multiple data exchange steps between the various products. With Aphelion Dev, generating a report can be done by only specifying the data that are needed (i.e., measurements and corresponding image). Then Dev automatically launches Microsoft Excel®, runs the default template or user macro, and generates the analysis report. Alternatively, the Report Generation task can either export grid values to a text or Excel file, or the image along with the associated overlay to the Windows® clipboard or in a file.



Upgrade Path for Advanced Capabilities

Aphelion™ Dev provides a comprehensive set of image functions that make possible the development of sophisticated image-based applications. Should your imaging project needs become more demanding, the Aphelion product family includes a large set of optional extensions that complement, enhance, and extend Dev's capabilities. Find out all the available extensions at <https://www.adcis.net/en/aphelion-extensions/>.

Product Distribution

Aphelion Dev is compatible with the most recent Microsoft® Windows® operating systems⁽²⁾, and cameras compliant with Microsoft DirectShow and IIDC 1394 Specification⁽³⁾. A version of its libraries, called Aphelion SDK, is also available for deployment including a basic image display capability.

Licensing

- Single User license (locked to a single computer system)
- Site license
- Multiple computer, non-concurrent use license (optional USB license key required)

Standard package

DVD containing the Aphelion Dev installation program, all extensions, acquisition device drivers⁽³⁾, multimedia demonstrations, and a User Guide hardcopy.

The Aphelion Dev installation program and user guide are also downloadable from <https://www.adcis.net>

⁽¹⁾ Processing and display of 3D images require two optional modules.

⁽²⁾ Refer to <https://www.adcis.net/en/aphelion-dev/> for the list of Windows version compatibility

⁽³⁾ Microsoft DirectShow protocol and IIDC 1394 Digital Camera Specification Version 1.31 as published by the 1394 Trace Association. Optional interfaces for digital cameras are also available (Refer to <https://www.adcis.net/en/applications/image-acquisition-device-interfaces/>)

⁽⁴⁾ Only available for BasicScript®



ADCIS S.A.
3, rue Martin Luther King
14280 Saint-Contest
France
Phone: +33 2 31 06 23 00
Fax: +33 2 31 06 23 09
www.adcis.net

Aphelion is a trademark of ADCIS S.A.
Microsoft, Windows, Excel, Visual Basic, Visual C#, and .NET are trademarks of the Microsoft Corporation. All rights reserved.

All other trademarks, product names, and company names used herein are for identification purposes only and may be trademarks or registered trademarks of their respective owners.